

vivaspaTM

Living the Healthy Life!



ADJUSTABLE INTENSITY
TO MASSAGE ACHING FEET



AIR COMPRESSION
TARGETS ANKLE TENSION



IMPROVES
BLOOD CIRCULATION



relax at
the push of
a button



kneading & heating shiatsu
360° AIR PRESSURE FOOT MASSAGER

MS-1600
USER'S MANUAL

Important Safety Precautions

When using your Kneading & Heating Shiatsu 360° Air Pressure Foot Massager, basic safety precautions should always be followed, including the following:

1. Read all instructions before using your device.
2. This appliance is not intended for use by persons (including children) with reduced physical sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
3. Stop use if the device causes you physical discomfort or abnormal feeling and consult a doctor.
4. Do not stand while using the massager.
5. Consult a doctor before use if you are pregnant, have a pacemaker, or any other health concerns.
6. Do not operate any appliance with a damaged cord or plug or after the appliance malfunction or has been damaged in any manner. Return appliance to manufacturer for examination, repair, or adjustment.
7. Do not use this device if you are inebriated or have taken medication that impairs judgment.
8. Do not let cord hang over the edge of table or counter, or touch hot surfaces.
9. Do not place on or near a hot gas or electric burner, or in a heated oven.
10. The heat function can cause the unit to become hot. Exercise caution when using the heat function and allow the unit to completely cool between uses.
11. Do not cover the unit with blankets or other coverings when in use.
12. Do not use outdoors.
13. Do not use the device where aerosol products are being used or where oxygen is being administered.
14. Keep out of reach of children.
15. Please recycle or dispose of your appliance properly based on the laws and rules of your municipality. Contact local recycling facilities and/or the manufacturer of your appliance for further information.

Introduction

Thank you for purchasing the Kneading & Heating Shiatsu 360° Air Pressure Foot Massager, item MS-1600. Relax after a long day with an at-home massage. The air compression targets ankle tension, the adjustable intensity massages your feet and improves circulation.

Package Contents

- MS-1600 Kneading & Heating Shiatsu 360° Air Pressure Foot Massager
- User's Manual and 1 Year Warranty Card

Features

- Wireless Remote Control
- Washable Foot Sleeves
- LCD Display & Massage Nodes
- Portable & Stylish Dome Shaped Design

Specifications

Rated Power: 40W

Rated Voltage: 110-250V

Timer: 5-30 minutes

A Closer Look



Powering On

1. Plug the power cord into a suitable outlet.
2. When the power is connected, a Red Indicator Light will illuminate.
3. Unplug the power cord when the massager is not in use.



Red Indicator Light



CAUTION

Do not squeeze, stretch, or damage the power cord.

Using Your Foot Massager

Control Panel



The Timer Indicator shows how much time is left on the massage. You can use the Timer Button just below it to pick how long you want the massage to last.

The Mode button allows you to toggle between three different massage modes.

The Intensity button allows you to toggle between different intensities ranging from low to high.

The Heat button turns on or off the heating function. The light will illuminate red when the button is 'on.'

The Power button turns the massager on or off.

Use the Wireless Remote to access all the above settings at your convenience.

Using Your Massager (Continued)

Read all instructions carefully before beginning use:

1. Sit down in a relaxed manner with the massager placed in front of you.
2. Press the Power button to turn on the massager. It will automatically be set for a 15-minute massage.
3. Place feet in the massager and relax.
4. Toggle between the different settings to your preference. It is not recommended to exceed a 30-minute massage.



Switching Modes

By pressing the Mode Button on the massager or on the remote you can choose among three different massage modes. Each mode offers varying treatment depending on your needs.

MODE 1


This mode provides the most complete massage experience. You will feel the massage pads rolling, kneading and gently scratching your heels, front soles and arches. In addition you will feel some air compression treatment for your ankles.

MODE 2

This mode focuses primarily on treatment for the arches and front soles with only light treatment for the heels.

MODE 3

This mode consists primarily of treatment for the ankles with heavy emphasis on compression.

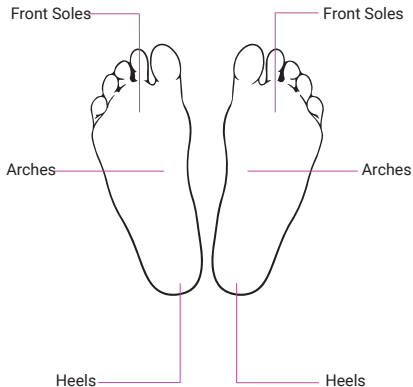

Maintenance & Care

After use, while the appliance is off, wipe with a soft cloth.

Keep the appliance away from extreme temperatures.

When not in use, keep the appliance in a safe place away from children or pet's reach.



CAUTION

DO NOT expose this appliance to extreme temperatures.

Do not let children play with this appliance.

Always Follow All Manufacturer's Instructions When Using.

Troubleshooting

Problem	Possible Causes	Possible Solutions
Unable to Boot	Power is off.	Turn on Power button.
	15-minute timer ran out.	Turn off and then turn on.

FCC Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Vivitar One Year Warranty

This warranty covers the original consumer purchaser only and is not transferable.

This warranty covers products that fail to function properly UNDER NORMAL USAGE, due to defects in material or workmanship. Your product will be repaired or replaced at no charge for parts or labor for a period of one year.

What Is Not Covered by Warranty

Damages or malfunctions not resulting from defects in material or workmanship and damages or malfunctions from other than normal use, including but limited to, repair by unauthorized parties, tampering, modification or accident.

To Obtain Warranty Service and Troubleshooting Information:
Call 1-800-592-9541 or visit our website at www.vivitar.com.

To receive Warranty service along with the name and address of an authorized product service center, the original consumer purchaser must contact us for problem determination and service procedures. Proof of purchase in the form of a bill of sale or receipted invoice, evidencing that the product is within the applicable Warranty period(s), MUST be presented in order to obtain the requested service. It is your responsibility to properly package and send any defective products along with a dated copy of proof of purchase, a written explanation of the problem, and a valid return address to the authorized service center at your expense. Do not include any other items or accessories with the defective product. Any products received by the authorized service center that are not covered by warranty will be returned unrepared.

Vivitar One Year Warranty (continued)

NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OTHER THAN THOSE EXPRESSLY DESCRIBED ABOVE SHALL APPLY. DISTRIBUTOR FURTHER DISCLAIMS ALL WARRANTIES AFTER THE EXPRESS WARRANTY PERIOD STATED ABOVE. NO OTHER EXPRESS WARRANTY OR GUARANTY GIVEN BY ANY PERSON, FIRM OR ENTITY WITH RESPECT TO THE PRODUCT SHALL BE BINDING ON DISTRIBUTOR. REPAIR, REPLACEMENT, OR REFUND OF THE ORIGINAL PURCHASE PRICE - AT DISTRIBUTOR'S SOLE DISCRETION -ARE THE EXCLUSIVE REMEDIES OF THE CONSUMER. IN NO EVENT WILL DISTRIBUTOR, ITS MANUFACTURERS, OR SAKAR INT. BE LIABLE FOR ANY INCIDENTAL, DIRECT, INDIRECT, SPECIAL, PUNITIVE OR CONSEQUENTIAL DAMAGES (SUCH AS, BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, BUSINESS, SAVINGS, DATA OR RECORDS) CAUSED BY THE USE, MISUSE OR INABILITY TO USE THE PRODUCT. EXCEPT AS STATED HEREIN, NO OTHER WARRANTIES SHALL APPLY. NOTWITHSTANDING THE FOREGOING, CONSUMER'S RECOVERY AGAINST DISTRIBUTOR SHALL NOT EXCEED THE PURCHASE PRICE OF THE PRODUCT SOLD BY DISTRIBUTOR. THIS WARRANTY SHALL NOT EXTEND TO ANYONE OTHER THAN THE ORIGINAL CONSUMER WHO PURCHASED THE PRODUCT AND IS NOT TRANSFERABLE.

Some countries, states or provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on warranties, so limitation or exclusions may not apply to you. This warranty gives you specific legal rights, and you may have other rights that vary from state to state or province to province. Contact your authorized retailer to determine if another warranty applies. This product is manufactured, distributed and sold by SAKAR International, Inc. All other trademarks are the property of the respective owner, who has not sponsored, endorsed or approved this product.

VIVITAR®