

INSTRUCTION
MANUAL

la parfait
COSMETICS

la parfait
L COSMETICS

INSTRUCTION MANUAL FOR THE **ULTRASCULPT MINI**
Skin care solutions custom tailored to you.

TABLE OF CONTENTS


ULTRASCULPT MINI DEVICE AND DISPLAY	1
DIRECTIONS AND KEY FEATURES	2
ULTRASCULPT MINI DEVICE MODES	3
WARNINGS STORAGE AND SPECIFICATIONS	4

800 689 4243

730 w. Hallandale Beach Blvd.
Suite #105, Hallandale Beach Fl,
33431


la parfait
COSMETICS

ULTRASCULPT MINI DEVICE AND DISPLAY

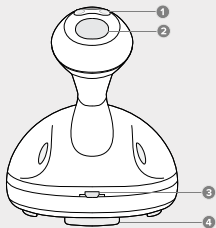

SIDE VIEW OF DEVICE

1	Power and setting buttons
2	Smart display screen
3	Charging port
4	Energy massage head

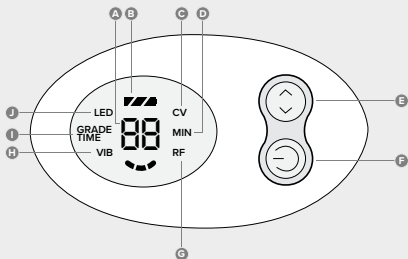

LCD SCREEN AND BUTTONS

A	Minute countdown
B	Battery life
C	CV treatment indicator
D	Minutes
E	Mode Intensity button
F	Power button / Cycle modes
G	RF treatment indicator
H	Vibration indicator
I	Treatment time / Intensity
J	LED treatment indicator

DIRECTIONS AND KEY FEATURES

DIRECTIONS

1. After bathing, if applicable, apply any of the recommended creams and/or Oils to the target area you are going to focus on.
 2. Press and hold the power button to turn on machine. (Default mode when the device is turned on is the Brighten skin mode). Press the power button to Cycle through the various modes and select the appropriate mode (see page 3).
- Note:** Look at the smart display to see what treatment indicators are showing up on the screen to make sure you are using the correct treatment!
3. Press the mode intensity button to cycle through the various intensities.
 4. Hold the handle firmly and begin massaging the skin in a circular motion. Working from the bottom to top.
 5. When your treatment has finished hold the power button to shut off the device.
 6. After using please clean the massage head with a soft lightly wet towelette.

KEY FEATURES

- **Ultrasonic High frequency vibration (CV)** Instantly wakes up the skin on a deep level improving and promoting cell metabolism and fat decomposition
- **Radio Frequency Shaping (RF)** Breaks up cellulite and takes advantage of the heating component to further shape, compact and sculpt the body
- **Wavelength 625nm high energy red light (LED)** Using the red LED and photosynthetic energy. You get the skin renewal and tightening effects of the Red LED while also promoting cell metabolism!
- **Mechanical Vibration (VIB)** An added layer of vibration helps accelerate the rate of fat decomposition. By contracting the collagen fibers.

ULTRASCLPT MINI

DEVICE MODES

DEVICE MODES

BRIGHTEN SKIN MODE

Uses LED Red Light and RF technology. The LED light activates the skin cells. Promoting cell metabolism, evening skin tones, and tightening the skin. A perfect environment for the RF heat to maximize its fat burning and body shaping magic.

Note: This mode is the default mode when booting on the machine. It has three levels of intensity. Slow bursts, rapid succession bursts, and lastly an always on option. Adjust based on your sensitivity.

FAT BURNING MODE

This utilizes both LED Red Light and Ultrasonic cavitation to loosen stubborn fat pockets, improve skin color and elasticity. Both LED and CV promote cell metabolism giving your skin the perfect push in the right direction! This mode only has one level of intensity.

HIGH BURST FAT MODE

Uses LED Red Light, RF technology, and mechanical Vibration. Using all three of these methods you take advantage of the cell metabolism and skin tightening from the LED lights in tandem with the Mechanical vibrations and RF heat gives the best environment to sculpt away those stubborn deposits

Note: This mode has three levels of intensity. Fast short bursts of light and vibrations, Slow bursts of LED with long intermittent vibrations, and lastly constant LED lights with long intermittent mechanical vibrations. Adjust based on your sensitivity.

SHAPING MODE

Utilizes both mechanical vibration and ultrasonic CV. Takes advantage of both technologies to loosen fat deposits and sculpt and tighten skin.

Note: This mode has three levels of intensity. Intermittent mechanical short vibrations. Intermittent mechanical long vibrations and mechanical long vibrations.

- Each treatment mode has a max time limit of fifteen minutes
-

WARNINGS

STORAGE AND SPECIFICATIONS

OPERATING WARNINGS

- Don't use if pregnant or have a skin disease
- Avoid using after a meal
- Don't exceed more than three times a week
- Don't overuse. Over using product will cause overheating and damage to the sensitive components
- Stop us and ask a doctor if any skin issues or irritation appear.
- Recommended for use on shaved skin.

STORAGE AND CARE

- Use with supplied charger. Don't use while plugged in.
- Not waterproof. Don't use in bath or sauna.
- Store in a dry place and keep out of reach from children.
- Clean and disinfect machine before and after each use.

Product	UltraSculpt Mini
Purpose / Intended use	Purpose / Intended use
Vibration frequency	680000 time/s
Unit size	90x90x110mm
Operating voltage	8.4V
Operating Current	0.14-2A
Charging voltage	5.0V
Power	16.8W

LAPARFAITCOSMETICS.COM