

INSTRUCTION
ULTRA SMOOTH IPL DEVICE

la parfait
COSMETICS



Instruction manual for la parfait's Ultra Smooth IPL Device
skin care solution custom tailored to you.


TABLE OF CONTENTS

WHAT'S INSIDE	1
DEVICE AND TREATMENT INSTRUCTIONS	2
SKIN AND HAIR COLOR TREATMENT GRAPH	3
POWER AND FREQUENCY TREATMENT CHART	4
USEFUL INFORMATION AND TIPS	5-6
WARNINGS AND INFORMATION	7
BEFORE AND AFTER PHOTOS	8

800 689 4243
2801 Greene street Bay 3A,
Hollywood, FL
33020

la parfait
COSMETICS


WHAT'S INSIDE THE BOX

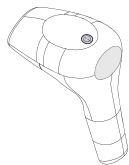

Ultra smooth IPL Device




Ultra Smooth Cooling Gel

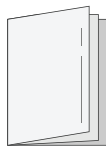

Instructional manual



Razor

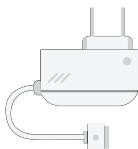

Charger










IPL Safety glasses

IPL DEVICE TREATMENT INSTRUCTIONS

- 1 Before anything else, shave hair from the skin surface that you want to treat.
- 2 Clean your skin thoroughly from any hair residue and apply the Ultra smooth cooling gel on the skin surface prior to use. (Can be also used without the gel on a clean dry skin)
- 3 Before you start to use the device, clean the flash window (quartz lamp)
- 4 Plug in the power adapter to a power source and plug the other end to your device. When connected the screen will show “HR”. Press and hold the power button for 3 seconds (the cooling fan will turn on and the display will show the remaining number of pulses).
- 5 A short press on the power button once at a time will switch between intensity level 1 to 8 and back down to 1 (it is recommended to follow the chart provided in this manual on page 4 for intensity levels for each area).
- 6 The device uses a proximity sensor built-in the flash, and the pulse will work only with skin contact. (Cover the flash window completely)
- 7 On the display make sure the “HR” is flickering, that means that the device is ready for the pulse. (If “HR” is not flashing the device will not emit light)
- 8 Press the Flash button, the device will flash one time. Lift device from the skin, wait till “HR” blinks again (takes 1 - 3 seconds) and repeat on area that was not treated. Alternatively a long press on the flash button will place the device into “auto mode” flashing automatically.
- 9 When treatment is finished, press the power button for 3 seconds to turn off the device.
- 10 Clean the flash window (quartz lamp) with an alcohol pad or moist wipes. (Do not use water on the device)

SKIN TONE

HAIR COLOR TREATMENT GRAPH

SKIN TONE	White	Ivory	Beige	Light Brown	Medium Brown	Dark Brown or Black
HAIR COLOR						
Black 	✓	✓	✓	✓	✓	✗
Dark Brown 	✓	✓	✓	✓	✓	✗
Brown 	✓	✓	✓	✓	✓	✗
Bark 	✓	✓	✓	✓	✓	✗
Blonde 	✗	✗	✗	✗	✗	✗
Red 	✗	✗	✗	✗	✗	✗
White 	✗	✗	✗	✗	✗	✗

Check your skin tone to find out if you can safely use the jpl ultra smooth device. Match the skin in the area you wish to treat to the range of the skin tones shown on the chart above. If your skin tone matches any of the skin tones labeled "not safe" for treatment, then do not use the device.

"✗" means not safe
"✓" means safe

BODY PART & SKIN TONE

LEVEL AND FREQUENCY OF USE CHART

BODY PARTS / SKIN TONE	WHITE	YELLOW	LIGHT BROWN	MEDIUM BROWN	BROWN
Armpit	(Levels 6-8) 3-4 times a week	(Levels 6-8) 3-4 times a week	(Levels 5-7) 2-3 times a week	(Levels 4-6) 2-3 times a week	(Levels 1-4) 2-3 times a week
Arm	(Levels 6-8) 3-4 times a week	(Levels 6-8) 3-4 times a week	(Levels 5-7) 2-3 times a week	(Levels 1-4) 2-3 times a week	(Levels 1-4) 2-3 times a week
Leg	(Levels 6-8) 3-4 times a week	(Levels 6-8) 3-4 times a week	(Levels 5-7) 2-3 times a week	(Levels 1-4) 2-3 times a week	(Levels 1-4) 2-3 times a week
Mouth Area	(Levels 4-8) 3-4 times a week	(Levels 4-8) 3-4 times a week	(Levels 4-6) 2-3 times a week	(Levels 1-4) 2-3 times a week	(Levels 1-3) 2-3 times a week
Face	(Levels 5-8) 3-4 times a week	(Levels 5-8) 3-4 times a week	(Levels 4-6) 2-3 times a week	(Levels 1-4) 2-3 times a week	(Levels 1-3) 2-3 times a week
Neck	(Levels 5-8) 3-4 times a week	(Levels 5-8) 3-4 times a week	(Levels 4-6) 2-3 times a week	(Levels 1-4) 2-3 times a week	(Levels 1-4) 2-3 times a week
Chest	(Levels 5-8) 3-4 times a week	(Levels 5-8) 3-4 times a week	(Levels 5-7) 2-3 times a week	(Levels 4-6) 2-3 times a week	(Levels 1-5) 2-3 times a week
Back	(Levels 6-8) 3-4 times a week	(Levels 6-8) 3-4 times a week	(Levels 5-7) 2-3 times a week	(Levels 4-6) 2-3 times a week	(Levels 1-5) 2-3 times a week
Bikini Line	(Levels 4-8) 3-4 times a week	(Levels 4-8) 3-4 times a week	(Levels 4-7) 3-4 times a week	(Levels 3-5) 2-3 times a week	(Levels 1-3) 2-3 times a week

IPL DEVICE USEFUL INFORMATION

USEFUL INFORMATION

- When using the device for the first time, it is recommended to check for any skin reaction by testing on the skin 2 hours prior to treatment. Test a small area under the upper arm. Start with a low level and increase gradually . Apply twice on the same area. If after 2 hours the area is still red, burned or in pain. Please consult your doctor before using the device again.
- If a slight discomfort occurs, use product again on a lower setting, however if discomfort continues consult your doctor before further use.
- After treatment, cool the treated area with ice or with a cold wet towel.
- It is strongly suggested not to use the device multiple times on the same part of the skin (up to 2 times for the same area).
- The device is suitable for all genders.
- The use of the IPL device is painless.
- It is safe for use on the face below the cheekbones, upper lip, chin, forehead and jawline.
- Avoid operating the device close or directly at the eyes.
- The hair growth cycle varies between individuals and could take up to 12-24 months, depending on the area that you treated. Since only hairs in their growing phase can be treated with light.
- It is recommended to do multiple treatments to achieve continuously smooth skin.
- We recommend following the treatment chart so you can achieve optimum results.
- Please avoid exposure to UV light sources, such as tanning booths or the sun.
- Avoid hot tubs, saunas and any skin care products that contains fragrance. Bleaching Creams and exfoliators, as they may cause irritation, such as itching and redness on the treated areas.
- While using the device, protect your eyes by using provided glasses.
- After each treatment, we recommend using our Body Revive Lotion.
- IPL is a highly effective treatment in killing the infection that causes fungal nails on both hands and feet by using Intense Pulsed Light (IPL). To achieve optimal results, we recommend a series of 3–4 treatments, spaced 1 month apart. It is essential to eradicate all fungal spores that may bloom after the first treatment. Improvement can usually be seen after one treatment.
- While treating fungal, keep using anti-fungal cream after each treatment. (keep using cream as directed)

IPL DEVICE INFORMATION AND BEST RESULTS

USEFUL INFORMATION

- Do not block the air vent, as the device will overheat and stop working.
- If device overheats, disconnect the power cord from the device and plug it back in.
- Hair that stays stuck on the device may cause blank spots but will not affect performance.
- As a result, due to the normal aging process of the tubes, the light bulbs will appear in different shades of black.
- Do not store the device near strong magnetic fields, damp environment or high temperatures.
- Handle with care and keep out of reach of children.
- Do not use any other light therapy treatments before and after using the IPL device.
- If you have any other issue or any questions please contact our customer service 800-689-4243

FOR BEST RESULTS

- On the first week we recommend using the device every other day (depends on the skin tone)
- On the second week, use twice on higher levels.
- After 4 weeks, use once a week at the highest level.
- If the hair keeps growing, continue to use more frequently.
- When the hair growth decreases, use the device every 2 months for maintenance.

THIS PRODUCT CAN BE USED FOR

- Excess facial hair grown on the upper lip, armpit, body, forehead, chest, back, legs, etc.
- Hair removal is suitable for dark hair tone such as black, dark brown, brown and light brown.
- After 4 weeks, use once a week at the highest level.
- It is not be suitable for white, gray, blonde, red or golden hair types.
- Has been proven to aid in fungal Treatment

WARNINGS AND INFORMATION

NOT SUITABLE FOR USE

- You have a known skin disorder.
- You have skin disease, such as allergic dermatitis, eczema, etc.
- You have photosensitization.
- If you are menstruating, pregnant, or lactating.
- Sunburn skin.
- Dark skin tone (Please refer to the skin tone hair color treatment graph on page 03)
- Broken skin, a scabbed or healed wound, cancer or hemangioma.
- Do not use on tattoos.
- Before and after treating the skin, avoid exposing treated area to the sun.

Lamp tube	Xenon quartz tube
Light source	Intense pulse light
Energy frequency	1.5-3.9j / cm ²
Wavelength	470nm-1100nm
Lamp tube life	900000 times
Adaptor	Ac 100V ~ 240V @ 12V DC 3A
Weight	210 grams
Size	168x112x78mm
Operating temperature range	5 ~ 30 c
Operating humidity range	Relative 30% ~ 60%

IPL DEVICE BEFORE AND AFTER

BEFORE AND AFTER RESULTS OF USING THE IPL DEVICE FOR HAIR REMOVAL TREATMENT



BEFORE AND AFTER RESULTS OF USING THE IPL DEVICE FOR FUNGAL TREATMENT



LAPARFAITCOSMETICS.COM