

Skin Rejuvenation Beauty Device for Face and Neck  
(USB charged)

Item No.: BD2

# GLO24K

UNLEASH THE POWER OF GOLD



**Important Reminder:** 1. Please charge to full before use.  
2. Please recharge the device between applications.

# Contents

Page

①

Product specifications

①

Instructions of the 3 modes

②

Patent

③

Product features

④

How to use on neck

④

How to use on face

④

How to clean

⑤

Cautions


⑥

Q&A


⑦

Technical information

# Product specifications

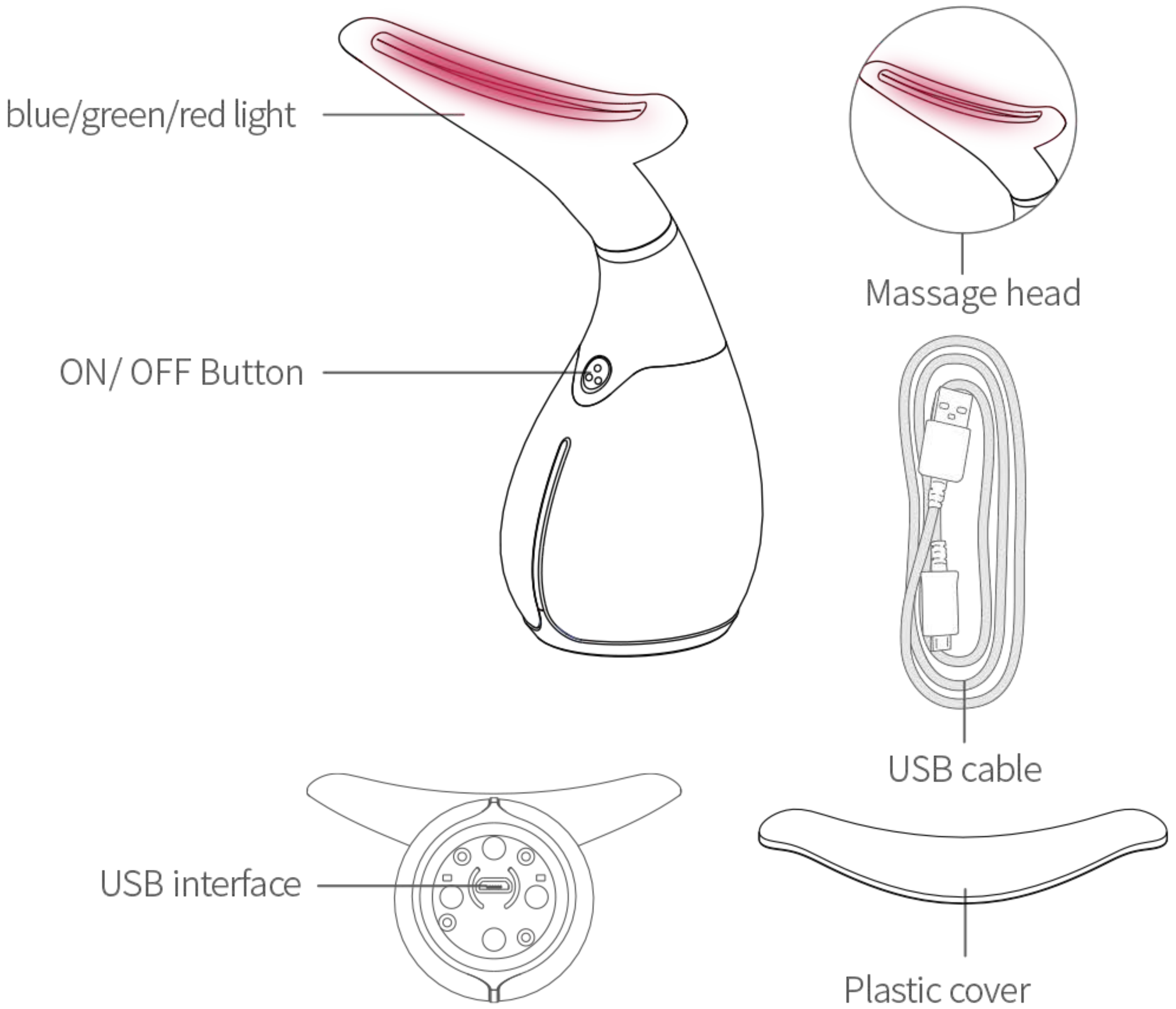


# Instruction of the 3 modes

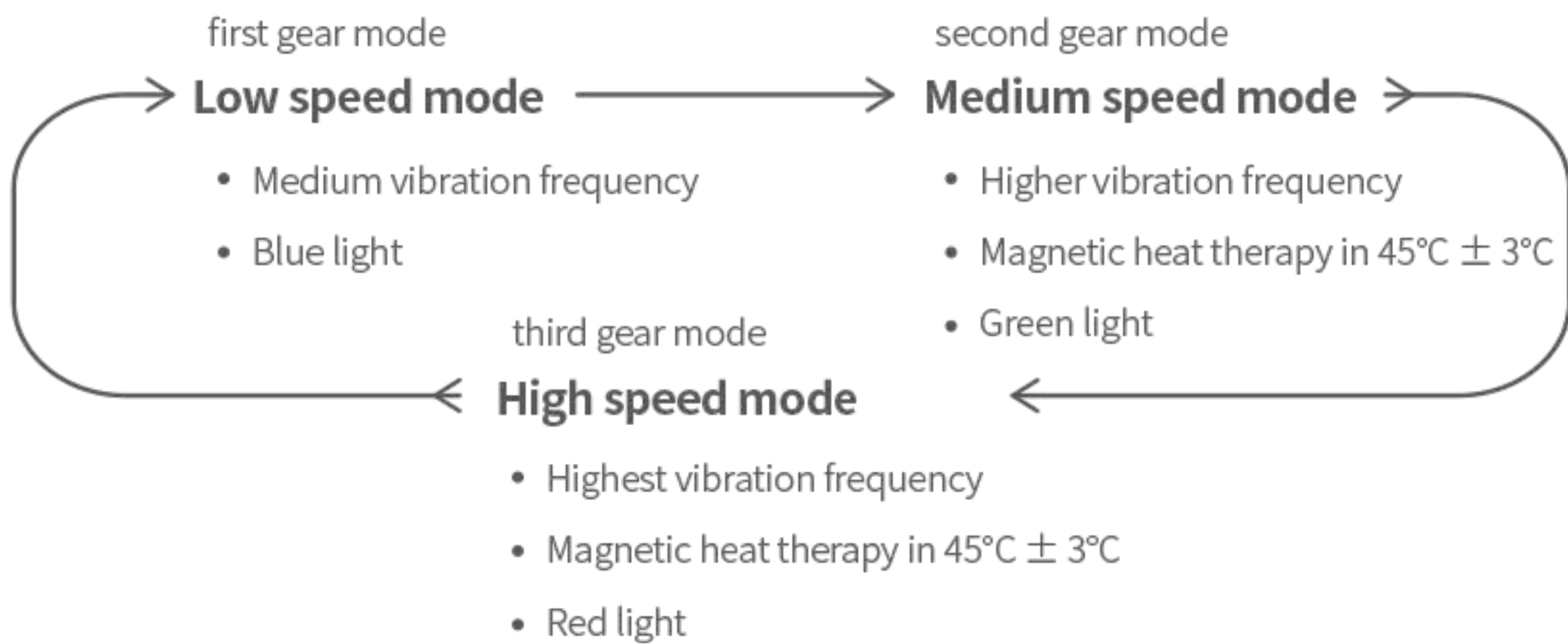



# Patent

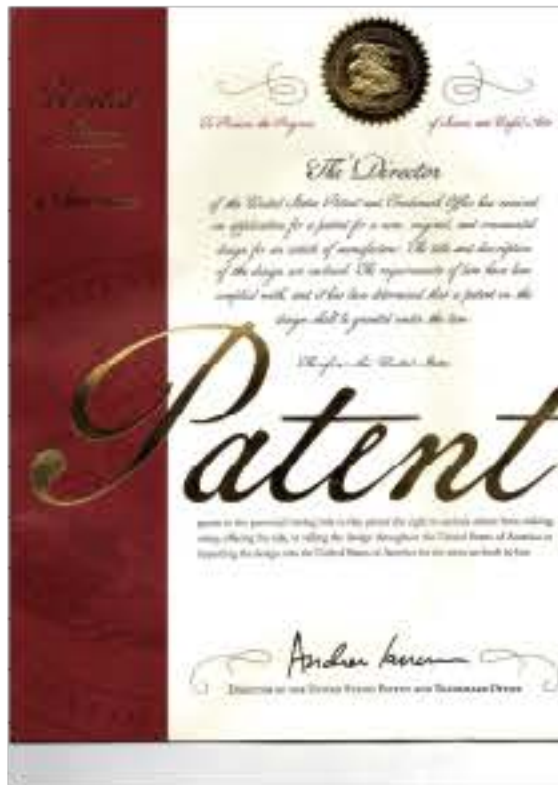

Appearance patent  
US D876,623 S



Appearance patent  
003539014-0001



Appearance patent  
30-0965355



Blue light is the coldest light, which can promote collagen synthesis, activate skin, improve sagging and refine pores. Also, it greatly contributes to improve oily and sensitive skin.



Green light can increase cell metabolism, promote the micro-circulation on skin surface and help with skin puffiness. Meanwhile, it has the effect of smoothing and balancing skin.





Red light can strengthen cell vitality, accelerate blood circulation as well as promote the growth of fibroblasts and collagen. It also has a significant influence on improving fine lines, wrinkles, sagging skin, flaccid, melanin and scars.



## Product features

- The massage head is specially designed for shoulder, neck and other parts of body, streamline outlook fit for ergonomic theory, can help transporting the energy to every part of skin.
- The massage heads heating to 45°C to activate the skin, matching with three IPL acoustic mode in low, medium and high speeds, and blue, green, red 3 color light massage modes, can promote the blood circulation, stimulate the neck acupuncture points, clear the meridians, improve circulation and promote metabolism, promote the discharge of long- accumulated toxic substances & waste, slow skin cell aging, to dilute the neck profile, less swelling, relax neck, make acupuncture treatments, rapidly eliminate fatigue, restore the role of physical fitness.
- Three sonic modes as low speed (blue), medium speed (green), high speed (red), choose mode according to personal needs.
- With mode memory function and 60 seconds intimate reminder function, it can easily adjust the massage time.
- USB cable charging, safety and environmental protection.
- Mini & portable, can be used anywhere any time.
- Newly adopted the bionics design concepts, shaped in dolphin tail, looks novelty and classic.

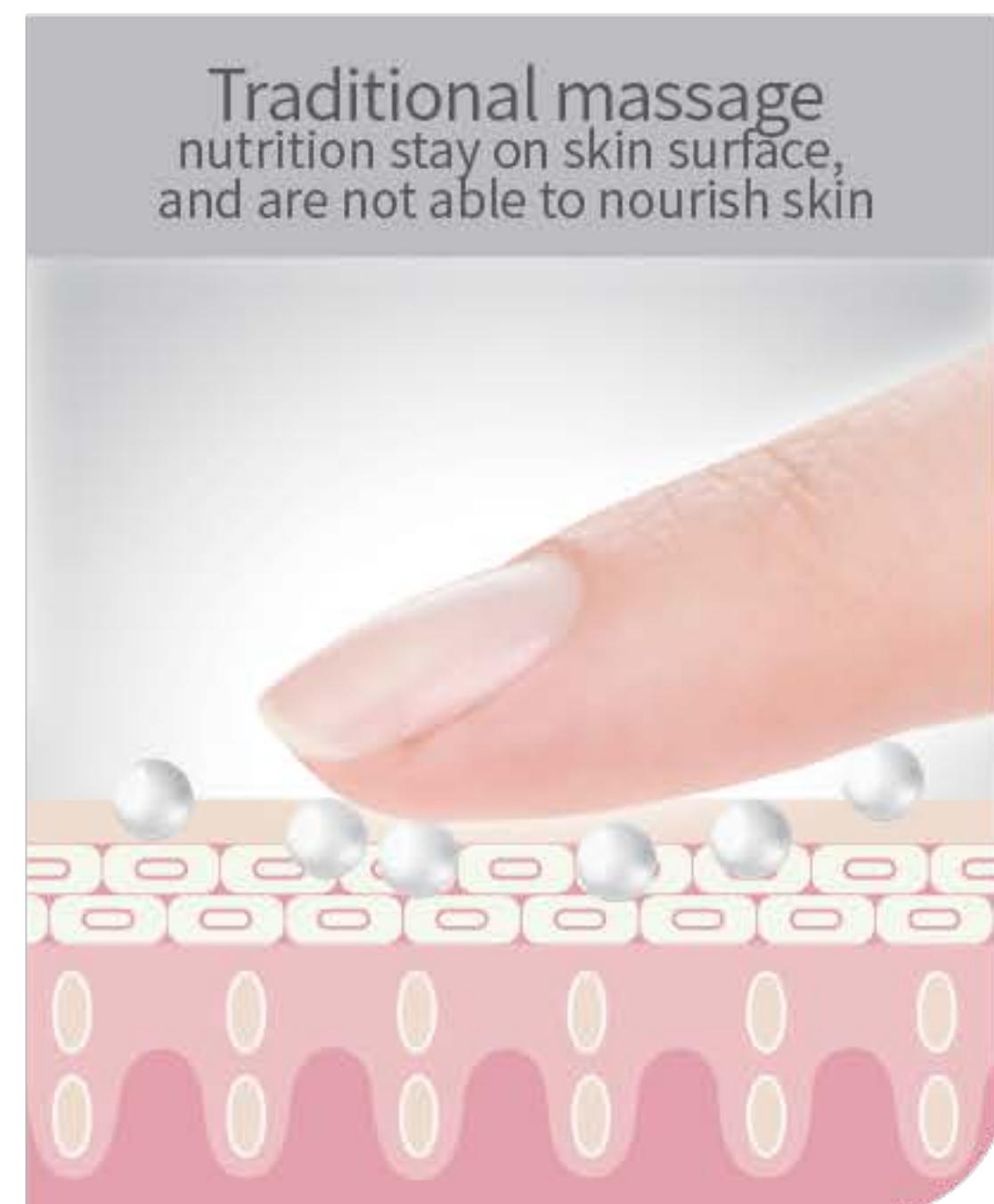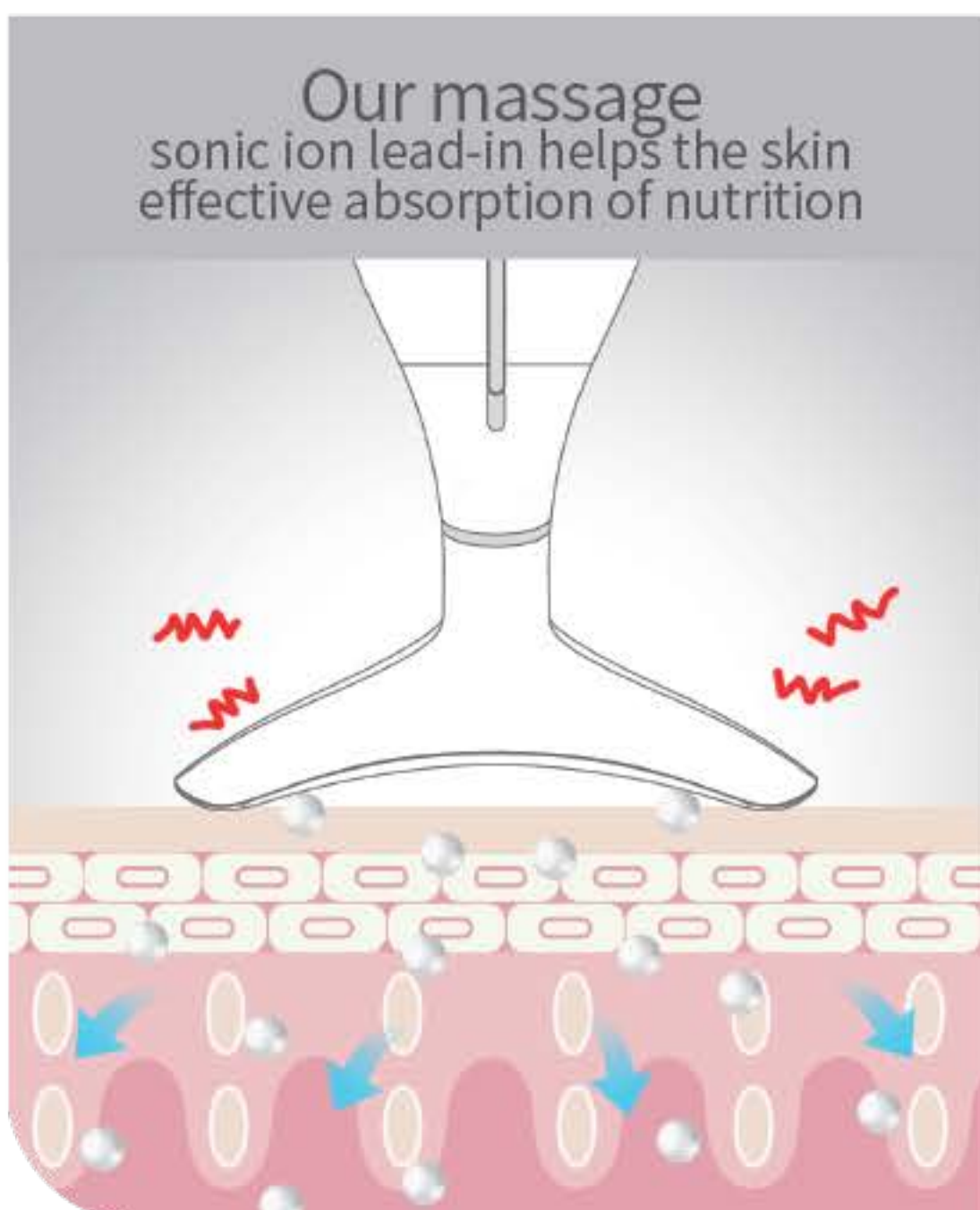




## How to use on neck

---

1. Apply a small amount of skin care cream on the neck area. Massage until fully absorbed.
2. Select your preference mode (blue/green/red) and apply the device on the skin in upwards/downwards motions.
3. Use the device on the back and front of the neck. Massage all acupuncture points in upwards/downwards motions. Use for 5-10 minutes at a time.

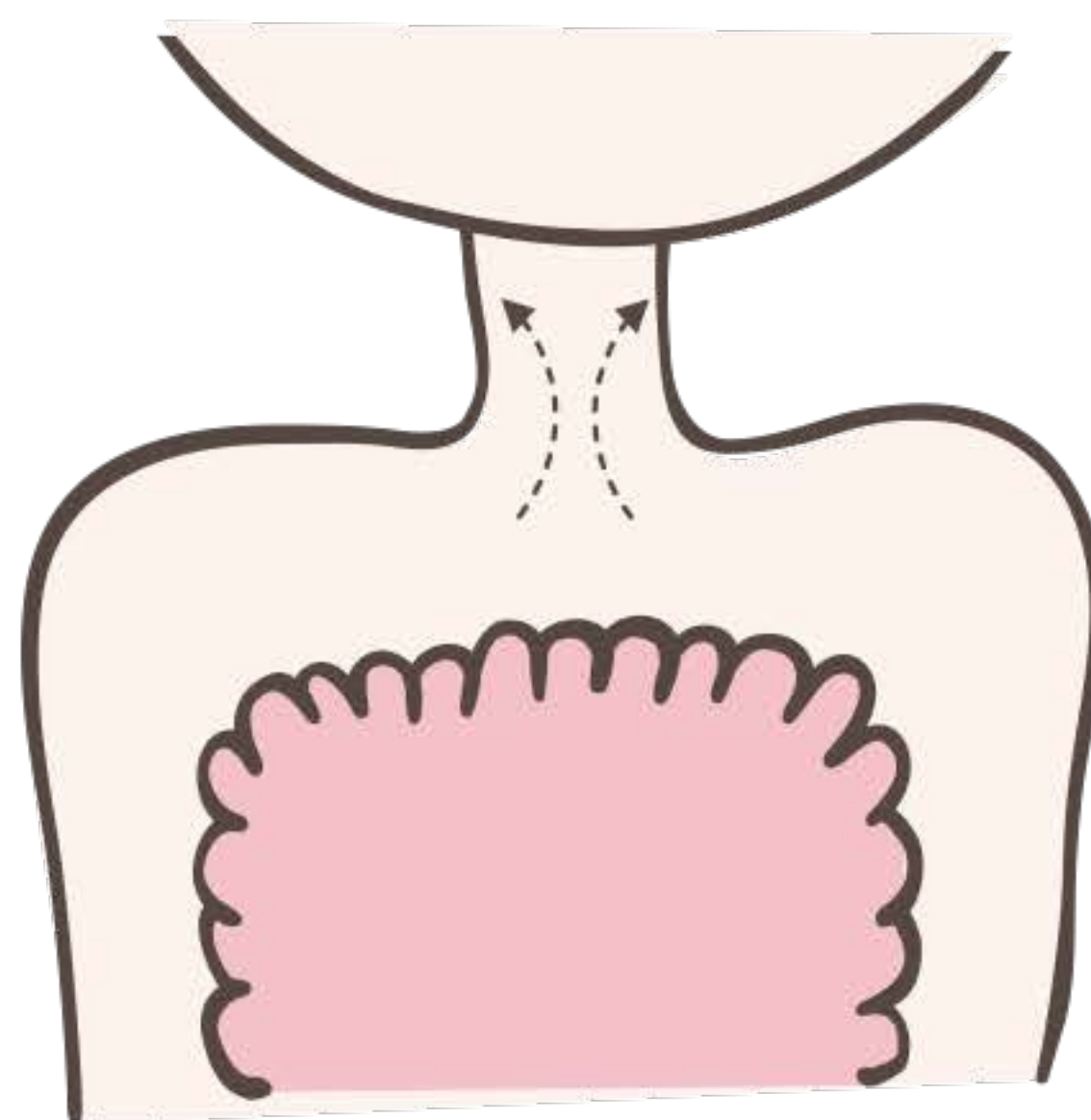

## How to use on face

---

1. Apply a small amount of skin care cream over clean skin. Massage until fully absorbed.
2. Select your preference mode (blue/green/red) and apply the device on the skin in upwards/downwards motions.
3. Massage all acupuncture points in upwards/downwards motions around the shape of your face. Use for 5-10 minutes at a time.



## How To clean

---


1. Press ON/OFF Button for one second to turn off the device.
2. Wipe the product with a dry cloth.
3. Close the plastic cover and place it in a ventilated dry environment.




## Cautions

---

1. Please read the product manual before using, and keep the manual for future use.
2. USB cable charging, charging from low battery station to saturation need about 2 hours, fully charged allows 2-12 hours of continuous use.
3. When battery is in shortage, flash in red light, when it is under battery charging, flash in blue light, when blue light is up on indicates charge in saturation.
4. The product has mode memory function, when switch it on, it will enter the mode you chose last time automatically.
5. It is not allowed to use this product when it is under charging.
6. The product has the function of heating, approximate temperature is in 45° C. Different people are in some different when sensing the temperature., some are likely to feel hot.
7. Please stop using the device if you feel uncomfortable or skin becomes irritated and contact the manufacturer.
8. The product is not waterproof; please do not place it in water.
9. Put the product on a dry and ventilated place, to avoid direct sunlight, and make sure no corrosive chemicals around.
10. The cleansing cloth using high-quality fleece material, can be washed repeatedly.
11. Friendly reminder: please use moderate skin care product. Excessive skin care product can not be thoroughly absorbed. In contrast, the leftover might possibly permeate into the device and damage it. Please clean the device in time with cleansing cloth after using.
12. After a long time of idle mode, the device will need about 15 minutes of charging to reactivate (indicator light will turn blue).
13. Keep out of reach of children, to avoid any danger.
14. Do not apply directly over eyes.

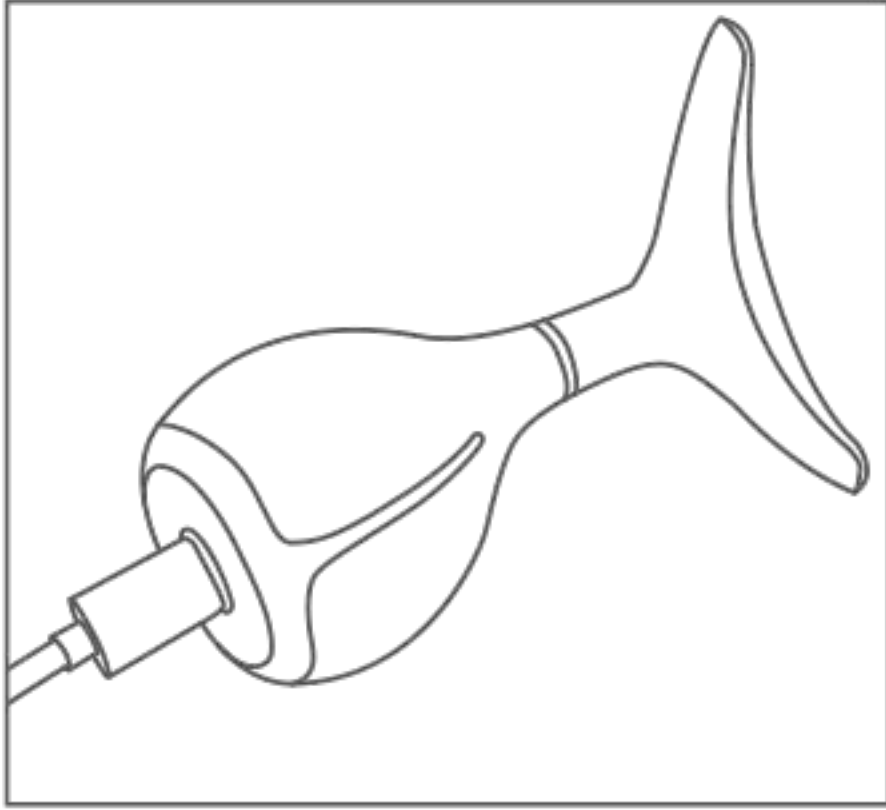

At the beginning of charging

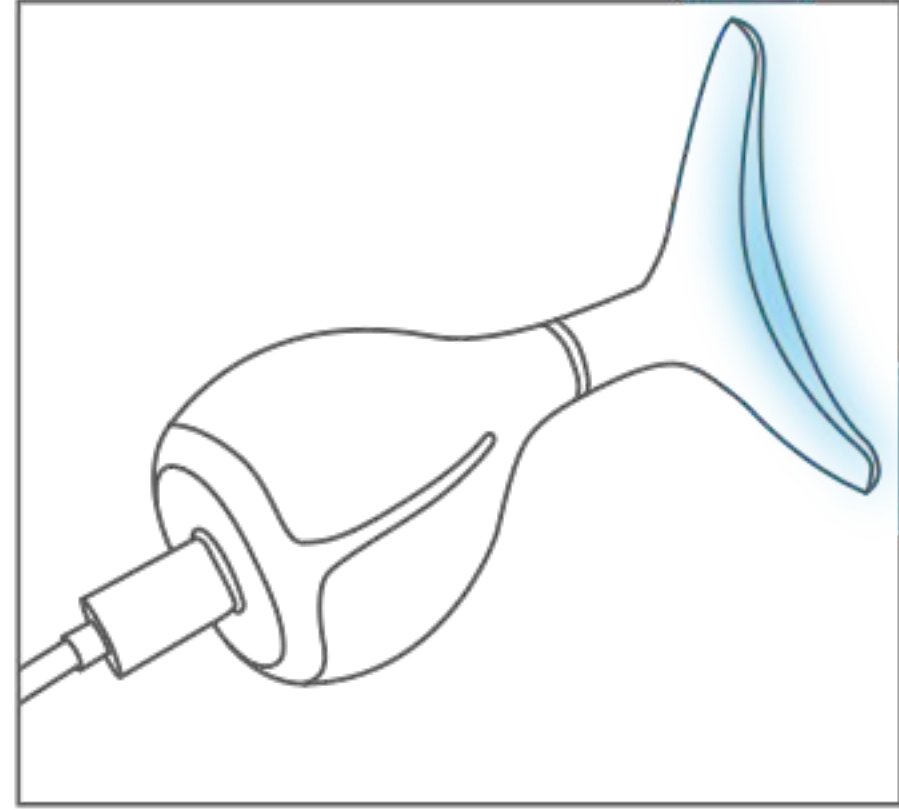

The product needs 15 mins to be the charging state

GLO24K offers a limited warranty on its beauty devices (check website for details).

## Q&A

### Question

The product suddenly has no vibration, does not work, what is the reason?

Can I use this product every day?

### Answer

- ▶ Check whether the product electricity is in shortage or not. Recharge device.

---

- ▶ Check if the device starts correctly.

---

- ▶ Restart the product after turn it off .

---

- ▶ Please contact the GLO24K to consult with the manufacturer.

---

- ▶ Yes, the vibration intensity and temperature adjustment of this device are intelligently designed by a professional team, can be used multi-times a day.



## Question

I suffer from allergies, can I use the device?

What is the recommended use of the device?

## Answer

Yes. This device is made with materials of the highest standards that are proven to be safe and not cause any skin irritation.

- ▶ However, in the rare case of skin irritation, allergy, or discomfort stop using the device and contact the manufacturer.

- ▶ The recommended use of the device is 5-10 minutes at a time any day you want.

## Technical information

---

Product name	Skin Rejuvenation Beauty Device for Face and Neck
Model	BD2
Rated Voltage	5V 
Rated Input Power	3W
Battery composition	Lithium
Battery type/ size	Lithium
Number of batteries	1
Battery weight ( grams )	19±0.3
Cells per battery	1
Watt hours per battery	2.4wh
Lithium content ( grams )	0.195

# GLO24K

UNLEASH THE POWER OF GOLD



**TO MAXIMIZE RESULTS SCAN HERE**

**For more information about the device please visit**

**[WWW.GLO24K.COM](http://WWW.GLO24K.COM)**

**Dist. by Y&E Ent., LLC**

**Atlanta, GA 30338 USA**

**71-75 Shelton Street, London, UK**